

seca 725

Mechanical baby scale with sliding weights



- Classic design and high-quality precision mechanics.
- Precise ascertainment of weight due to zero adjustment.
- Softly curved weighing tray for maximum safety.
- Easy-care and mobile.

seca 725:

The sturdy classic.

A baby gains weight quickly – its weight has almost doubled after four to six months. As soon as the baby begins to crawl and run, baby fat will quickly disappear. In order for the baby's development to be monitored safely from the very beginning, instruments are required in pediatrics in which one can rely on long-term. A baby scale that will work in all surroundings: The mechanical baby scale seca 725 has been around for over 50 years and is still being used by health professionals all over the world.



Classic design and high-quality precision mechanics.

The seca 725 combines the classic, flat construction with high-quality precision mechanics. With its smooth running weights, it is easy to use and provides precise results. The baby lies safe and secure in the softly curved weighing tray. The robust steel construction and cast levers ensure years of durability and accuracy. All this has made this mechanical baby scale a favorite for many decades.

Precise weight ascertainment due to zero adjustment.


The zero adjustment ensures the reliable determination of the net weight – a soft mat can simply be tared away, showing the actual weight of the baby. The baby scale has lbs and kg measurements on the balance beam and the weight graduations are very small, 1/4 oz. and the weight capacity is 32 lbs. These graduations are hard to beat on any mechanical baby scale you can find.

Easy-care and mobile.

The weighing tray is easy to remove from the base for cleaning. The hard-wearing, scratch-proof finish and the smooth, seamless surface makes the seca 725 extremely easy to clean.



To simplify cleaning, weighing tray can be removed.

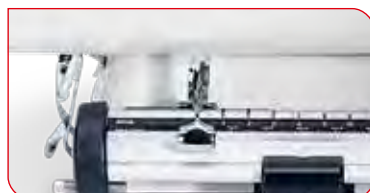

The reset-to-zero function guarantees reliable determination of net weight by simply deducting the extra weight of padding or diaper.

seca 725

Technical Data

- Capacity: 32 lbs, 16 kg
- Graduation: 1/4 oz, 5 g
- Dimensions (WxHxD):
21.7 x 7.1 x 11.4" /
550 x 180 x 290 mm
- Dimensions/tray (WxHxD):
21.7 x 4.3 x 10.6" /
550 x 110 x 270 mm
- Weight: 14.1 lbs / 6.4 kg
- Functions/Accessories:
reset-to-zero function